


I/O TTL INTERFACE FOR SMART CENTER (BN-SMART-IOCBL)

Access 8 Digital TTL inputs on the Smart Center

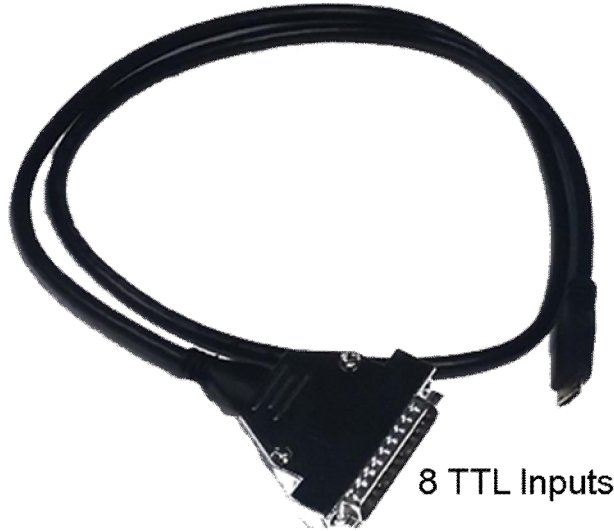

This 1 meter cable connects to the small I/O port on the Smart Center and terminates in a male DSub25 connector. Use to connect to E-Prime, SuperLab, or other visual presentation systems.

Pin assignments match E-Prime style pin assignments:

- the first TTL digital input channel is pin 2 of the DSub
- the final TTL digital input is pin 9 of the Dsub

Note: The Smart Center does not presently support digital output, however this functionality may be added in a future firmware and software revision. Contact BIOPAC to indicate your interest in Smart Center digital output lines.