


SS1LA SHIELDED ELECTRODE ADAPTER

The fully-shielded electrode interface cable permits high resolution recording of biopotential signals. The 3-meter adapter cable accepts standard 1.5 mm female Touchproof connectors. Use this lead adapter with:

- LEAD120 and EL120 Contact Post Electrodes
- LEAD110 Series and
 - EL160 Series Reusable Gold Cup Electrodes
 - EL250 Series Reusable Ag-AgCl Electrodes
 - EL350 Series Bar Electrodes
 - EL450 Series Needle Electrodes
- LEAD140 Series Clip Leads

To use disposable electrodes, connect 2 x LEAD110S-R/W shielded electrode leads to the VIN+ and VIN- inputs, and 1 x LEAD110 unshielded lead to the ground input. The LEAD110 series electrode leads will interface with any BIOPAC disposable electrode.

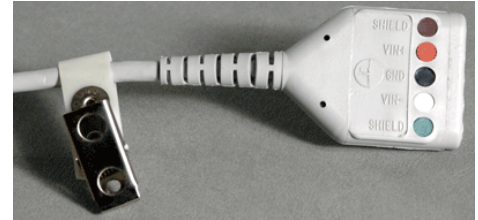

SS1LA SPECIFICATIONS

Cable length	3-meter
Termination	standard 1.5 mm female Touchproof connectors

Note: The SS1L is a 3-meter electrode adapter for older style 2 mm pin connections. To convert 2 mm pin connections to Touchproof 1.5 mm connections, use CBL201.