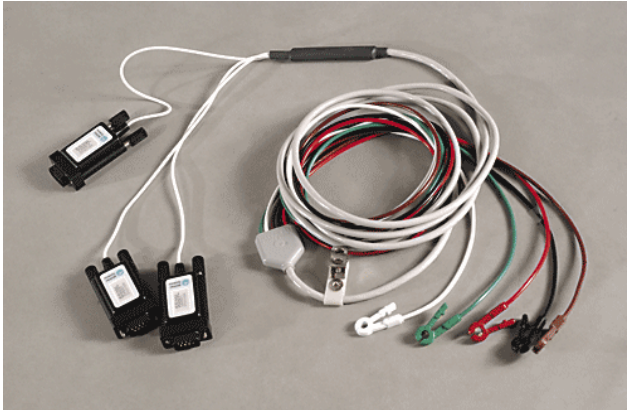


SS29L MULTI-LEAD ECG CABLE

The SS29L Multi-Lead ECG Cable permits high-resolution ECG recordings. This multi-lead set can simultaneously record Leads I, II, III, aVR, aVL, aVF, plus one precordial chest lead V(1-6). A 12-Lead ECG recording can be obtained by alternating the chest lead electrode from position V1 through V6. The cable terminates in three Smart Sensors that connect to the MP3X.

SS29L SPECIFICATIONS

Input Cable Length:	2 meters
Electrode Lead Length:	1 meter
Internal connection:	Built-in Wilson terminal
Electrode interface:	Connects to standard snap-connector disposable electrodes (EL503)